




## BOVIE<sup>®</sup> A952-V ELECTROSURGICAL GENERATOR

The Bovie A952, is a High Frequency Desiccator with cutting ability. You might call it a split personality. One half of the unit is a high frequency desiccator, an upgrade to our very popular Aaron<sup>®</sup> 940. The other half is a Bovie Generator, giving you 50 watts of cutting power and coag capabilities – plenty of power for most cutting procedures performed in the veterinarian's office.

- High Frequency Desiccator with 50 watts of Cutting capabilities
- Two machines in one convenient package
- Adjustable in 1/10th watt increments in the first 10 watts, for the most accurate output available
- No grounding pad required when working in the desiccation mode

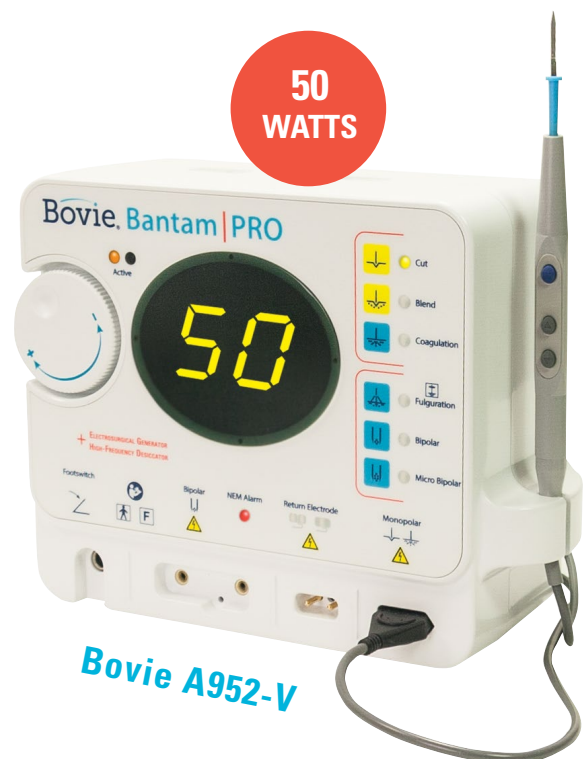

# BOVIE® A952-V, ELECTROSURGICAL GENERATOR



## STANDARD EQUIPMENT

- A952 (1) electrocautery generator,
- A902 (1) handpiece
- A910-P (10 pk) non-sterile disposable handpiece sheaths\*
- A910ST (2) sterile disposable handpiece sheaths\*
- ES01R (1), ES20R (1) reusable electrodes
- A833 (1) reusable electrodes
- A1254C (2) reusable return electrode cord
- A1204P (1) grounding plate
- A837 wall mount kit\*
- 110 VAC Hospital-grade power cord\*
- User's guide (MC-55-238-001)\*

\* Not Pictured



A902



ES01R, ES20R



A833



A1254C



A1204P

## COMPLEMENTARY PRODUCTS †



A813  
Table Top Stand †



A803  
Footswitch †



A827V  
Forcep Cable †



A824  
Reusable Bipolar Forceps †

† Sold Separately



Reusable Loop  
Electrodes †



Reusable Electrodes †



Tungsten Wire  
Loop Electrodes †

## PRODUCT SPECIFICATIONS

### DIMENSIONS

Height: 188mm (7.4 in.)  
Width: 228mm (8.98 in.)  
Depth: 105mm (4.13 in.)  
Weight: < 2.26 kg (< 5 lbs.)

### OUTPUT POWER

Cut: 50 watts  
Blend: 50 watts  
Coagulation: 50 watts  
Fulguration: 40 watts  
Bipolar: 40 watts  
Micro Bipolar: 40 watts

### OUTPUT FREQUENCY

304 – 451kHz

### LINE FREQUENCY

50Hz

### LINE VOLTAGE

A952: 100-240 VAC

### WARRANTY

Four Years