

Bovie®

SPECIALIST|PRO

Electrosurgical Generator



The Bovie® Specialist Pro is the electrosurgery solution for in demand physician offices and specialty facilities. Designed to operate like an operating room generator in a compact lower watt package. The multipurpose unit is ideal for monopolar and bipolar office surgical procedures.

The Specialist|PRO offers 5 common energy outputs controlled by the new sleek, push-button design, with up to 120 watts of power. The modalities include Cut, Blend, Coagulation, Fulguration and Standard Bipolar.

Efficiency, Versatility, and Safety ...
Designed benefits for you and your patients

Industry-standard connections allow for use of a variety of monopolar and bipolar accessories providing energy versatility across a variety of specialties such as OBGYN, Urology, Plastic Surgery, Dermatology, ER and more.

Bovie NEM™ neutral electrode monitoring increases patient safety by ensuring proper pad connectivity during procedures.

The Bovie Specialist|PRO is backed by the unmatched **4-year manufacturer warranty**, offered only by Bovie.



Visit BovieMedical.com/specialistpro for more information on the Bovie Specialist|PRO.

BOVIE® SPECIALIST|PRO ELECTROSURGICAL GENERATOR



STANDARD EQUIPMENT

- A1250S electrosurgical generator,
- ESP1 sterile disposable pencil
- ES01 blade, ES02 needle, ES20 ball (1 each)
- sterile disposable electrodes
- ESRE (5) split disposable return electrode
- A1252C (2) reusable grounding cords
- 110 VAC Hospital-grade power cord*
- 220 VAC Hospital-grade power cord*
- User's guide (MC-55-237-001)*

* Not Pictured



ESP1



ES01, ES02, ES20



ESRE



A1252C

COMPLEMENTARY PRODUCTS †



SE02
Smoke Evacuator †



BV-1253B
Monopolar Footswitch †




A1255A
Footswitch Adapter †




Disposable LLETZ/LEEP
Loop and Ball Electrodes †



ESREC - Split Disposable
Return Electrode with cord †

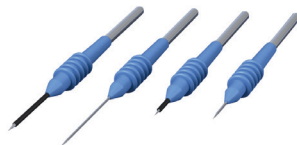

Disposable SuperCut™
Needles Electrodes †

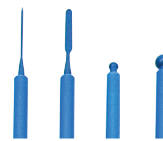

Disposable Resistick II™
Electrodes †



Disposable Dermal Loops †

† Sold Separately

PRODUCT SPECIFICATIONS

DIMENSIONS

Height: 15.2cm (6 in.)
Width: 26cm (10.25 in.)
Depth: 30.5cm (12 in.)
Weight: < 4Kg (< 9 lbs.)

OUTPUT POWER

Cut: 120 watts
Blend: 90 watts (60% Cut, 40% Coag)
Coagulation: 80 watts
Fulguration: 40 watts
Bipolar: 30 watts

OUTPUT FREQUENCY

350 – 600 kHz

LINE VOLTAGE

A1250S: 100-240 VAC

LINE FREQUENCY

50 – 60 Hz

WARRANTY

Four Years