


Sentinel Series

Scout Series



Color-Coding Options - Ease of Inventory Management

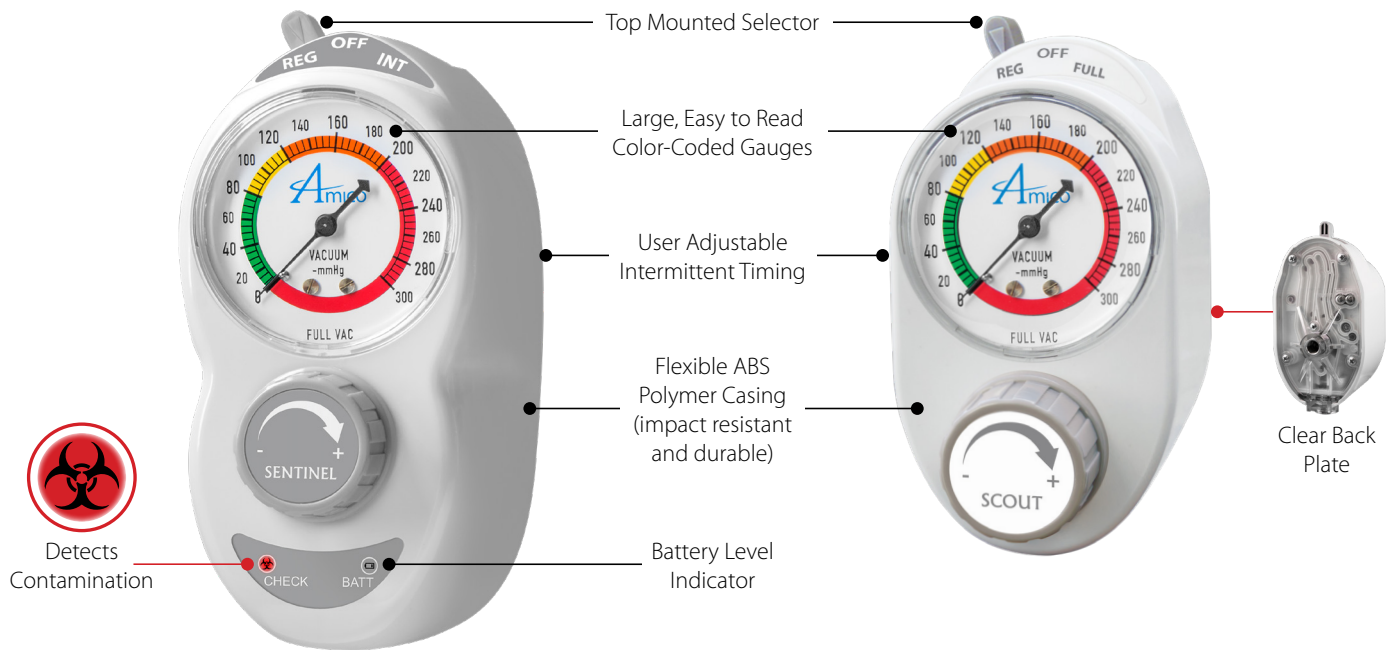


By Department:

- ED
- Pediatrics
- ICU

By Model:

- Surgical/Continuous High
- Continuous/Intermittent
- Pediatrics
- Neonatal

Continuous Models - White Labels

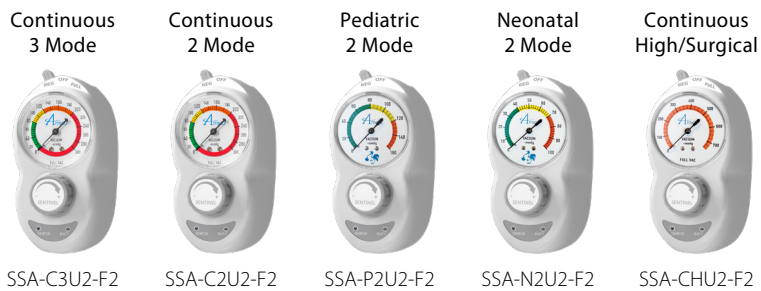

Continuous/Intermittent Models - Grey Labels



Features

- (Sentinel Only) Only regulator that detects contamination in **REAL TIME**
- (Scout Only) Only regulator on the market with a disposable, see-through back body
- Disposable body eliminates need for cleaning
- Helps to eliminate the dangers of over-suctioning by virtue of design
- Quick and easy maintenance - serviced in under 2 minutes
- Full bodied color-coding options for inventory management
- Ability to upgrade from Analog to Digital gauge
- Ability to reverse flush
- Warranty Periods:

Territory	Sentinel	Scout
North America	12 Years (Analog & Digital)	10 Years (Analog) 12 Years (Digital)
International	5 Years (Analog & Digital)	5 Years (Analog & Digital)

How Contamination Spreads

- **37% of regulators** in hospitals are **contaminated***
- Contaminants include “growth of well-established nosocomial pathogens” **known to cause Hospital Acquired Infections (HAI)***
- Contaminants **spread from a regulator to the wall-side canister and back to a patient’s stomach within 24 hours***

Identify and Remove Contamination

Step 1 Contamination is Detected with **Sensor** or **Clear Back Plate**



Step 2 Remove Contaminated Body




Step 3

Reverse Flush with a Cold Sterilant or Discard the Contaminated Body



Step 4 Reassemble with only 1 NEW Component

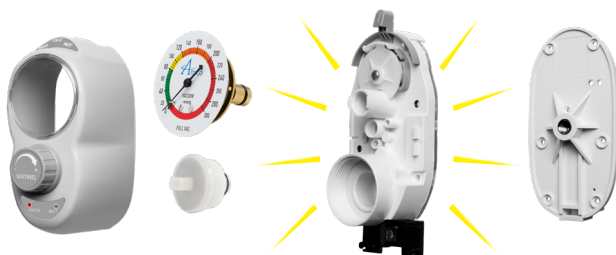

Step 5 Reset the Regulator (Sentinel Only)



Roll Stand for Vacuum Regulators



- Overall height 30" (76.2 cm)
- 5 swivel casters – 2" (5.08 cm)
- Stainless steel body
- Top handle for easy mobility
- Side hook for excess tubing or hoses
- 2 mounting locations for suction canister placement
- Regulator mount (VAC DISS Male) – available for single or dual regulators

Part#: S-RS-30-DISS-S*

*Canisters, holders, tubing, and suction regulator not included

Accessories



Vacuum/Safety Traps




Disposable Canisters and Holders



Autoclavable Canisters and Holders


Model Numbers

SXY-ZZXY-ZZ(X)*

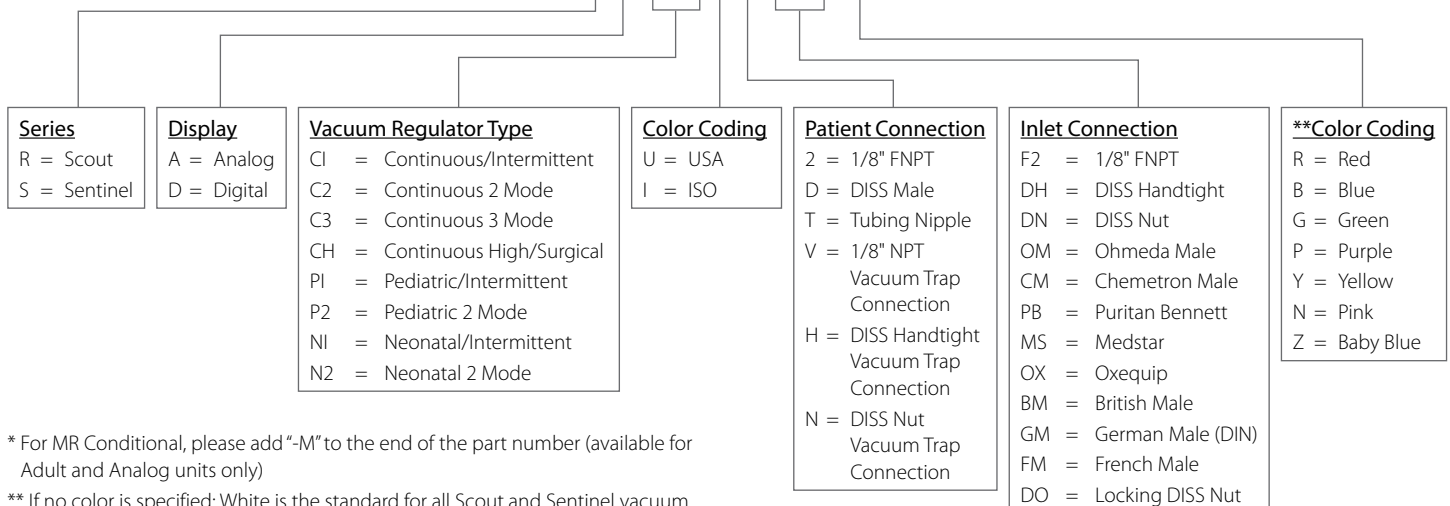

* For MR Conditional, please add "-M" to the end of the part number (available for Adult and Analog units only)

** If no color is specified: White is the standard for all Scout and Sentinel vacuum regulators

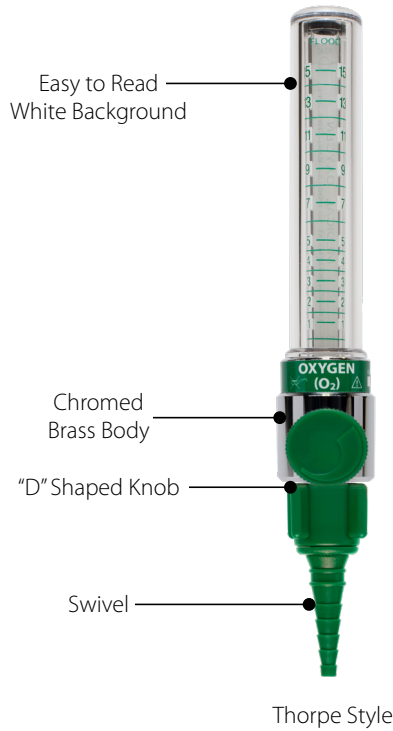

Dial Style



Dual/Y-block
Flowmeter Assembly

Features

- Complete line of Flowmeters from Neonatal to Adult
- Standard 1/8" FNPT connection, but also customizable with different adapters and accessories
- MRI options are available
- Available in USA and ISO colors
- Smallest and lightest chromed, brass bodied Flowmeter in the industry
- Color coded inner tube and labeling
- Easy to service
- No strip "D" shape knob
- Includes an inlet filter and inlet extension
- 5 year warranty North America (3 year International)

Accessories



Reusable Humidifier



Power Takeoff



Chain
(Flowmeter to Outlet)




Metal Swivel with Chain



Swivels/Christmas Trees

Model Numbers

FMX-XXY-XX(YY)*

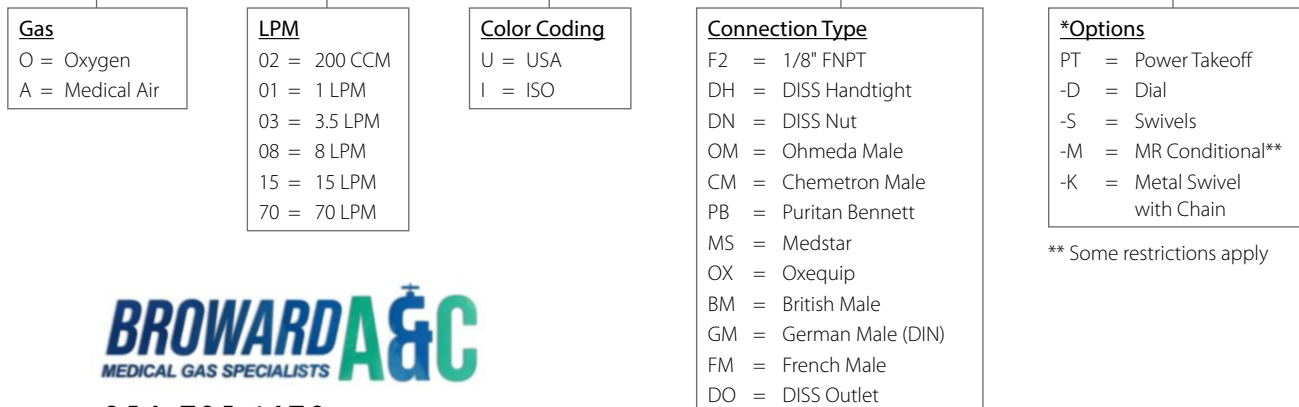

954-725-1470 Ext. 403

